

# Double Bed Raisers MF 2/4

Multifit Hospital Supplies Ltd - Enquiries Ph 0800 102 803 - Fax 0800 102 800  
E Mail : info@multifit.co.nz Website www.multifit.co.nz

**Unique design links 3 bed blocks together making the raised bed safer and more stable for your client.**



- **Tested** to British Standards / Stability of furniture ) to ensure stability of the raised bed and the safety of your client
- **Cost Effective** - durable construction means they can be reissued over and over.
- **Ease of installation** eliminates the need for you to get under the bed to fit individual blocks .
- **Rubber pad** improves slip resistance making the bed more stable.
- **No tools required** to clamp the blocks in position.
- **Powered Bed Jack** for trouble free lifting of any size bed is also available from Multifit

#### Method of installation

Just go to the foot of the bed, slide out the double bed raiser so that the blocks align with the legs of the bed. Clamp in this position by tightening the wing nuts on the underside.

Using a Bed Jack, lift the bed and slide the three connected blocks under the centre of the bed from the bed side. Because you have pre positioned the blocks, they will be automatically aligned with the bed legs without further adjustment. Press on the foot of the bed and the bed will rock upward, now fit the raiser at the head of the bed. Finally lift the foot of the bed and fit the Double Bed Raiser to the last three legs.

All this is accomplished without getting under the bed, putting your back at less risk when fitting the raiser and giving your client a safer raised bed .




*Double bed raiser gives a really solid base to the bed , no matter how high it is raised*

#### Specifications

Max User weight	380 kg
Length adjustment	1150 - 1700 mm
Sold in sets of	1,2 or 3
Non Slip foot pad size	90 x 90 mm
Lip around foot pad	20 mm high
Construction	Furniture grade pine and colour coated sheet metal
Std heights available	50,75,75 & 100mm
Finish	Clear Lacquer

# Multifit update

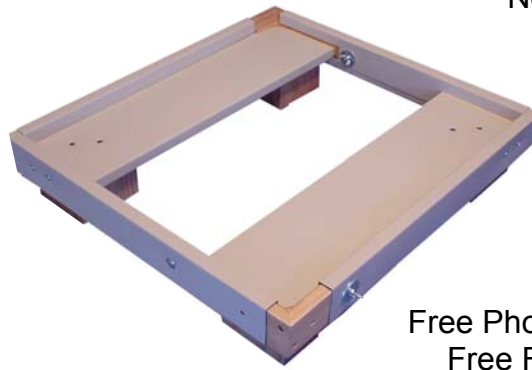

**Your recognized specialist in Chair Raisers**

Free Phone 0800 102 803  
Free Fax 0800 102 800  
E Mail : [info@multifit.co.nz](mailto:info@multifit.co.nz)  
Website : [www.multifit.co.nz](http://www.multifit.co.nz)

## OT workshops very worth while



**Beatrix van de Klundert**  
*OT from Dunedin reports on  
the OT Workshops held  
recently in Cromwell*

It was a very worthwhile opportunity and many thanks to Multifit for my sponsorship to the OT workshops in October.

There were many great workshops to attend and I can honestly say that all were worthwhile. There are some very innovative initiatives happening in our communities.

Some of the workshops I attended included Occupational Therapists taking up the

challenge of working in PHO's which proved a challenging topic as to the type of roles and areas we could work in. This included looking at mental health issues, older person's health, public health, education and identifying new initiatives to develop in our changing health climate.

As always with any occupational therapy workshop it is standard to participate. This included working in small groups, developing new ideas and always presenting back to the groups.

Occupational Therapists never seem to get an opportunity to sit back and be idle! It was a wonderful opportunity to meet new people and to put faces to names especially all the people from ENABLE funding, Accessable and the Ministry of Health who attended. This was a golden opportunity to attend workshops presented by Emily Gooder, Sue Primose and Rachel Tatham to discuss potential changes to Enable funding and Accessable especially in regard to the accredited assessor process. There were some challenging questions and an opportunity for therapists to put their views forward and to potentially change some of these processes to hopefully benefit the assessors and also the clients. All the workshops were well attended with some turned away if they had not enrolled for them previously!

Regards

*Beatrix van de Klundert*  
*Sponsored to the OT workshops by Multifit*  
*See [www.multifit.co.nz](http://www.multifit.co.nz) to find out what else Beatrix had to say about the workshops*

## Double & Single Bed Raisers are now on the Enable common list



### Double Bed Trough

Code MF 2/4  
Adjustment 1150 - 1700mm  
Heights 50,75,100 & 150mm  
Available in sets of 2 or 3

### Single Bed Trough

Code MF 2/2  
Adjustment 815 - 1100mm  
Heights 50,75,100 & 150mm  
Available in sets of 3 only

Enable have just finished their tender rounds for standard equipment and I am pleased to say that we now have two additional raisers on the Common list. Both the Single and Double Bed Raisers have been added as feedback from OT's suggested that they are stronger and safer than other methods of raising a single or double bed. Multifit products on the Enable Common list are :

Raiser	Height in mm	Code
Chair Raiser Standard	50,75,100&150	MF 1/1
Chair Raiser Lazy Boy	50,75,100&150	MF 1/3
Chair Raiser Big Softie	50,75,100&150	MF 1/7
Sofa Bed Trough (set of 2)	50,75,100&150	MF 2/1
Single Bed Trough (set of 3)	50,75,100&150	MF 2/2
Double Bed Trough set of 2 & 3	50,75,100&150	MF 2/4

## Hints & Tips

A lot of the couches today have big fat feet. We can make sofa raisers with the troughs wider. Simply take a measurement of the foot if it's square, or the diameter if it's round. Put a piece of paper under the foot and draw around it. You can then measure across the circle to get the diameter. Ask for a quote as it's not a common list

# JENZ SAY



The clocks are forward . We have summer in our sights. It's amazing how suddenly we feel energized when the sun is shining. A real shame that so many elderly people just don't make the effort to soak up a few "rays" It is so beneficial.

I had a great time in Australia. See me and "some Charlie" I found in Sydney! Fantastic place in the winter. Did a fair amount of travel. Sydney and up to Noosa. Got some great holiday snaps, caught up with family, played golf (very badly), drank too much wine, maxed out the Visa. Just the usual stuff!

Have caught up at work now so get those phones ringing.

What about a Xmas Raiser for the "old folks that have everything"?

We also have smaller items they could find useful in the car like the Handy Bar.

Look forward to your calls.

Cheers Jennifer

## Handy Bar


Code MF HB



**Makes a perfect gift**

### Additional hand hold makes transfers

**safer.** Simply insert the Handybar into the door latch U-shaped striker on the door frame of your vehicle. The U-shaped striker supports the Handybar while the handle extends out from your car providing you with an additional handhold.

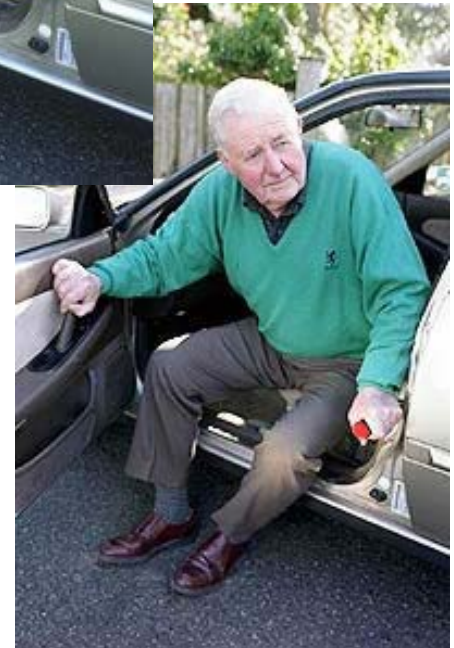

- Provides an additional handhold when getting in or out of your car, truck or van.
- Strong and stable - will give you added confidence and support making your transfers in and out of your vehicle safer.
- Fits most vehicles (both driver and passenger side-doors) see [www.handybar.com](http://www.handybar.com) for exceptions
- Increases independence and mobility
- Max user weight 158 kg—350 lbs
- Non Slip Handle
- Convenient size and weight

### Emergency Escape

As an added bonus the Handybar has a built in window breaker and seat belt cutter. This could protect you and your family from entrapment in a flooding or smoking vehicle. For more info : [www.handybar.com](http://www.handybar.com)

# Your easy order form

Free phone 0800 102 803

Free fax 0800 102 800

Multifit Hospital Supplies Ltd

Web site : [www.multifit.co.nz](http://www.multifit.co.nz)

E Mail : [info@multifit.co.nz](mailto:info@multifit.co.nz)

Name : ..... Date .....

Address .....

Phone ..... Fax .....

Delivery Address : ..... Assessor No : .....

Order  Trial  Quote  2-3 days delivery

Please invoice.....

Client ID	Order Code	Product Description	Height of unit	Quantity

## Win Morning Tea for your department

We really value your support. You promote our business in many ways and we appreciate it. To be in the monthly draw to win Morning Tea for your department tick the box and fax the form back to us.



## Last months' winners

**Aug** Annette Lane Takapuna Hospital  
**Sept** Anna Whiteside Balclutha

## Double handled bed rail

A really good solution for people who have minimal range of arm movement



Professionally tested for strength and stability

## Multifit Guarantee

- The equipment will retain its original function for **at least 2 years** or we will replace it **free of charge**.
- The equipment will fit or we will customise one that will fit.
- Delivery within **5 working days** or we will shout your **whole OT department morning tea**. ( customised & developmental equipment may take a little longer )